


# Scented Pre-filled Inflatable Disposable Anaesthesia Face Masks

Product Data Sheet

## PRODUCT: Scented Pre-filled Inflatable Disposable Anaesthesia Face Masks

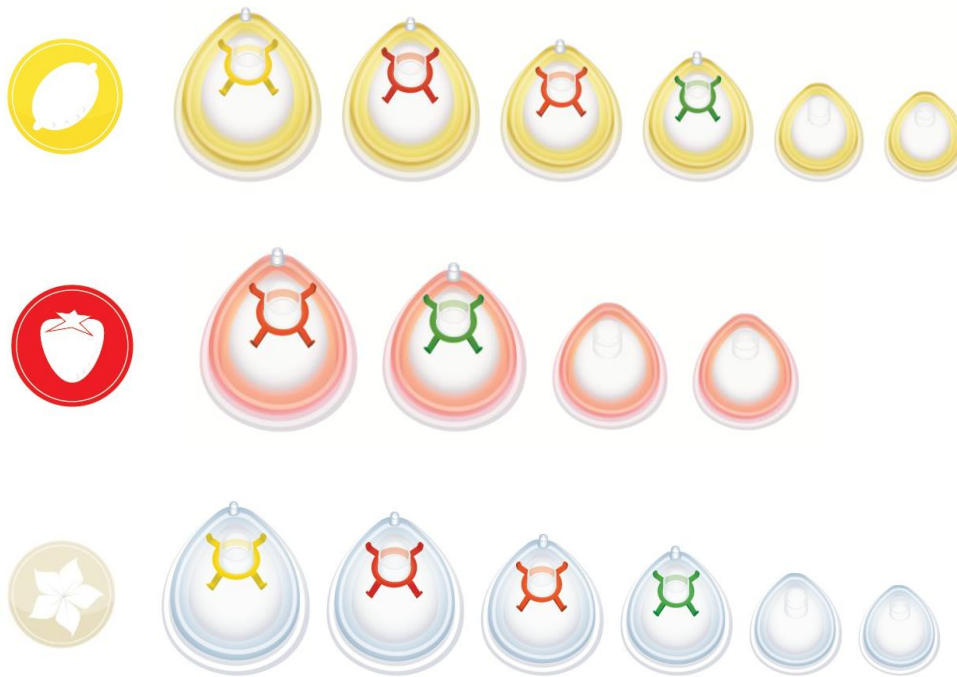

### DESCRIPTION / FUNCTION

Flexicare's Scented Pre-filled Inflatable Disposable Anaesthesia Face Masks are designed to channel anaesthesia gases or mixtures of oxygen to the patient for inhalation via the mouth and nose. Anaesthesia masks are to be used in conjunction with anaesthesia and ventilator breathing systems. They are scented as some patients receiving anaesthesia, particularly children, can find the experience unpleasant and a selection of three aromas can overcome this anxiety.





### WARNINGS



- Do not use if packaging is open or damaged
- Single use device only. Do not re-use, sterilise or autoclave
- After use, this product is contaminated waste and may be a potential biohazard. Handle and dispose of in accordance with hospital policy and applicable local guidelines and regulations.





### SPECIFICATION / KEY FEATURES

- Lemon, Strawberry or Vanilla scented
- Soft inflatable cushion can be inflated or deflated to enhance the seal to suit the patient's anatomical contours
- Clear face piece allows visual observation of patient status
- Compatible with all standard catheter mounts and breathing systems with 22mm & 15mm connectors
- Colour coded hooks to indicate mask size

## PRODUCT RANGE &amp; PART NUMBERS

Product Code	Description			Connector
038-51-400SFL	Lemon Scented Inflatable Disposable Face Mask	Large Adult	Size 5 	22mm
038-51-410SFL	Lemon Scented Inflatable Disposable Face Mask	Adult	Size 4 	22mm
038-52-415SFL	Lemon Scented Inflatable Disposable Face Mask	Large Child	Size 3 	22mm
038-52-420SFL	Lemon Scented Inflatable Disposable Face Mask	Medium Child	Size 2 	22mm
038-52-430SFL	Lemon Scented Inflatable Disposable Face Mask	Small Child	Size 1	15mm
038-53-440SFL	Lemon Scented Inflatable Disposable Face Mask	Neonatal	Size 0	15mm

038-52-415SFS	Strawberry Scented Inflatable Disposable Face Mask	Large Child	Size 3 	22mm
038-52-420SFS	Strawberry Scented Inflatable Disposable Face Mask	Medium Child	Size 2 	22mm
038-52-430SFS	Strawberry Scented Inflatable Disposable Face Mask	Small Child	Size 1	15mm
038-53-440SFS	Strawberry Scented Inflatable Disposable Face Mask	Neonatal	Size 0	15mm

038-51-400SFP	Vanilla Scented Inflatable Disposable Face Mask	Large Adult	Size 5 	22mm
038-51-410SFP	Vanilla Scented Inflatable Disposable Face Mask	Adult	Size 4 	22mm
038-52-415SFP	Vanilla Scented Inflatable Disposable Face Mask	Large Child	Size 3 	22mm
038-52-420SFP	Vanilla Scented Inflatable Disposable Face Mask	Medium Child	Size 2 	22mm
038-52-430SFP	Vanilla Scented Inflatable Disposable Face Mask	Small Child	Size 1	15mm
038-53-440SFP	Vanilla Scented Inflatable Disposable Face Mask	Neonatal	Size 0	15mm

## MATERIALS

Component	Material
Cushion	Polyvinyl Chloride (PVC)
Cushion Mask Shell	Polyvinyl Chloride (PVC)
Non-return Valve	Acrylonitrile Butadiene Styrene (ABS) / Silicone / Polyvinyl Chloride (PVC) / Polyoxymethylene (POM) / Stainless Steel
Hook	Polypropylene (PP)

## LATEX CONTENT

Flexicare Scented Anaesthesia Face Masks do not contain natural rubber Latex.

## DEHP CONTENT

Flexicare Scented Anaesthesia Face Masks do not contain phthalate DEHP.

## SINGLE USE

Flexicare Scented Anaesthesia Face Masks are single use devices and are to be disposed of after one use.  
DO NOT Reuse, Sterilise or Autoclave.

## STERILITY

Flexicare Scented Anaesthesia Face Masks are supplied non-sterile.

## STORAGE

Store in a cool, dry place out of direct sunlight.

## SHELF LIFE

Flexicare declares a shelf life of five years from the date of manufacture. This is based on the stability of the devices' components and raw materials sourced. Expiry date is clearly marked on individual product pouch.

## DISPOSAL CONSIDERATIONS

Dispose as clinical waste, in accordance with hospital policy, local guidelines and regulations.

## PACKAGING MATERIALS

Primary - Polybag Polyethylene (PE)  
Secondary - Carton Cardboard, Quantity per Carton = 50