

# Suction Catheters

Product Data Sheet

## PRODUCT: Suction Catheter



### DESCRIPTION / FUNCTION

Flexicare Suction Catheters are designed for removal of excess fluids from the oral, tracheal and endobronchial areas within a patient using vacuum suction.








Flexicare Suction Catheters include Standard and Vacuum Controlled. The range is essentially designed for suction implementation, and will connect to the suction source via Flexicare's Suction Tubing.

### WARNINGS














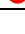
- Do not use if packaging is open or damaged.
- Single-use device only. Do not re-sterilise or autoclave.
- After use, this product is contaminated waste and may be a potential biohazard. Handle and dispose of in accordance with hospital policy and applicable local guidelines and regulations.

### SPECIFICATION / KEY FEATURES

- Smooth Tip at the patient end is rounded to reduce risk of injury
- Smooth Main Tube facilitates easy insertion through all types of airways, tubes and connectors
- Double Sided Eyes are smooth and free from sharp edges
- Clear Tube allows clinicians to visualise suctioning of patient
- Colour coded connectors for easy identification of tube diameter
- Supplied Sterile - Ethylene Oxide (EO) gas
- Conforms to requirements of ISO 8836 - Suction Catheters for use in the Respiratory Tract
- Catheter Tube Length 500mm
- Latex free
- DEHP free

Designated Size		Inside Diameter (mm)	Colour Identification
French (Charriere) Equivalent (F)	Outside Diameter (mm)		
6	2.00	1.00	Light Green 
8	2.66	1.75	Light Blue 
10	3.33	2.00	Black 
12	4.00	2.45	White 
14	4.66	2.95	Dark Green 
16	5.33	3.40	Orange 
18	6.00	3.90	Red 

**PRODUCT RANGE & PART NUMBERS**

Product Code	Description	Quantity
050-10-06	Suction Catheter without Vacuum Control, 6fr, 2.00mm 	100
050-10-08	Suction Catheter without Vacuum Control, 8fr, 2.66mm 	100
050-10-10	Suction Catheter without Vacuum Control, 10fr, 3.33mm 	100
050-10-12	Suction Catheter without Vacuum Control, 12fr, 4.00mm 	100
050-10-14	Suction Catheter without Vacuum Control, 14fr, 4.66mm 	100
050-10-16	Suction Catheter without Vacuum Control, 16fr, 5.33mm 	100
050-10-18	Suction Catheter without Vacuum Control, 18fr, 6.00mm 	100
050-20-06	Suction Catheter with Side Vacuum Control, 6fr, 2.00mm 	100
050-20-08	Suction Catheter with Side Vacuum Control, 8fr, 2.66mm 	100
050-20-10	Suction Catheter with Side Vacuum Control, 10fr, 3.33mm 	100
050-20-12	Suction Catheter with Side Vacuum Control, 12fr, 4.00mm 	100
050-20-14	Suction Catheter with Side Vacuum Control, 14fr, 4.66mm 	100
050-20-16	Suction Catheter with Side Vacuum Control, 16fr, 5.33mm 	100
050-20-18	Suction Catheter with Side Vacuum Control, 18fr, 6.00mm 	100

**MATERIALS**

Component	Material
Main Tube	Medical grade Polyvinyl Chloride (PVC)
Colour Coded Connector	Medical grade Polyvinyl Chloride (PVC)

## **LATEX CONTENT**

Flexicare Suction Catheters do not contain natural rubber Latex.

## **DEHP CONTENT**

Flexicare Suction Catheters do not contain phthalate DEHP.

## **SINGLE-USE**

Flexicare Suction Catheters are Single-use devices and to be disposed of after one use.

DO NOT Reuse, Re-sterilise or Autoclave.

## **STERILITY**

Flexicare Suction Catheters are supplied sterile by Ethylene Oxide (EO) Gas method.

## **STORAGE**

Store in a cool, dry place out of direct sunlight.

## **SHELF LIFE**

Flexicare declares a shelf life of five years from the date of sterilisation. This is based on the stability of the devices' components and raw materials sourced. Expiry date is clearly marked on individual product pouch.

## **DISPOSAL CONSIDERATIONS**

Dispose as clinical waste, in accordance with hospital policy, local guidelines and regulations.