

Pre-filled Inflatable Disposable Anaesthesia Face Masks

Product Data Sheet

PRODUCT: Pre-filled Inflatable Disposable Anaesthesia Face Masks



DESCRIPTION / FUNCTION

Flexicare's Pre-filled Inflatable Disposable Anaesthesia Face Masks are designed to channel anaesthesia gases or mixtures of oxygen to the patient for inhalation via the mouth and nose. Anaesthesia masks are to be used in conjunction with anaesthesia and ventilator breathing systems. They are compatible with all standard catheter mounts and breathing systems with 22mm & 15mm connectors.



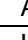
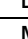
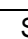
WARNINGS

- Do not use if packaging is open or damaged
- Single use device only. Do not re-use, sterilise or autoclave
- After use, this product is contaminated waste and may be a potential biohazard. Handle and dispose of in accordance with hospital policy and applicable local guidelines and regulations.

SPECIFICATION / KEY FEATURES

- Clear face piece allows visual observation of patient status
- Compatible with all standard catheter mounts and breathing systems with 22mm & 15mm connectors
- Soft inflatable cushion can be inflated or deflated to enhance the seal to suit the patient's anatomical contours
- Available in seven sizes, from neonatal to extra large adult
- Colour coded hooks to indicate mask size

PRODUCT RANGE & PART NUMBERS

Product Code	Description	Connector
038-51-390	Pre-filled Inflatable Anaesthetic Face Mask Extra Large Adult Size 6 	22mm
038-51-400	Pre-filled Inflatable Anaesthetic Face Mask Large Adult Size 5 	22mm
038-51-410	Pre-filled Inflatable Anaesthetic Face Mask Adult Size 4 	22mm
038-52-415	Pre-filled Inflatable Anaesthetic Face Mask Large Child Size 3 	22mm
038-52-420	Pre-filled Inflatable Anaesthetic Face Mask Medium Child Size 2 	22mm
038-52-430	Pre-filled Inflatable Anaesthetic Face Mask Small Child Size 1	15mm
038-53-440	Pre-filled Inflatable Anaesthetic Face Mask Neonatal Size 0	15mm

MATERIALS

Component	Material
Cushion	Polyvinyl Chloride (PVC)
Cushion Mask Shell	Polyvinyl Chloride (PVC)
Non-return Valve	Acrylonitrile Butadiene Styrene (ABS) / Silicone / Polyvinyl Chloride (PVC) / Polyoxymethylene (POM) / Stainless Steel
Hook	Polypropylene (PP)

LATEX CONTENT

Flexicare Pre-filled Inflatable Anaesthesia Face Masks do not contain natural rubber Latex.

DEHP CONTENT

Flexicare Pre-filled Inflatable Anaesthesia Face Masks do not contain phthalate DEHP.

SINGLE USE

Flexicare Pre-filled Inflatable Anaesthesia Face Masks are single use devices and are to be disposed of after one use.

DO NOT Reuse, Sterilise or Autoclave.

STERILITY

Flexicare Pre-filled Inflatable Anaesthesia Face Masks are supplied non-sterile.

STORAGE

Store in a cool, dry place out of direct sunlight.

SHELF LIFE

Flexicare declares a shelf life of five years from the date of manufacture. This is based on the stability of the devices' components and raw materials sourced. Expiry date is clearly marked on individual product pouch.

DISPOSAL CONSIDERATIONS

Dispose as clinical waste, in accordance with hospital policy, local guidelines and regulations.

PACKAGING MATERIALS

Primary - Polybag Polyethylene (PE)

Secondary - Carton, Quantity per Carton = 50