

Advanced
Solutions
in Acute
Respiratory
Care



Criterion[®] OxiCheck[™]


Portable,
Removable
Inline
O₂ Analyzer


The Criterion[®] OxiCheck[™] is a compact analyzer designed to “spot check” oxygen concentration.

Connected in seconds directly to the 22mm output connector of a Whisperflow[®] Generator or most ventilators, OxiCheck can be used to confirm the FiO₂ output being delivered to a patient quickly and accurately.

The OxiCheck is portable and can be carried in your pocket, ideal where a quick check of oxygen delivery is required.


Simple to use and calibrate, the OxiCheck also preserves battery life by automatically switching off after 90 seconds. The device lasts for 12 months and requires no maintenance.

- Low cost performance
- Weighs less than 100 grams
- Easy to calibrate
- Automatic power saving feature
- Battery included which operates for the life of the analyzer
-  Latex Free

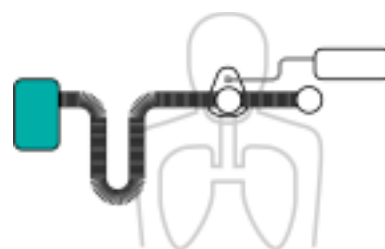

Criterion® OxiCheck™

OxiCheck™ is a compact, portable, inline O₂ analyzer. It is convenient to insert and remove from breathing circuits and works in conjunction with WhisperFlow Generators.

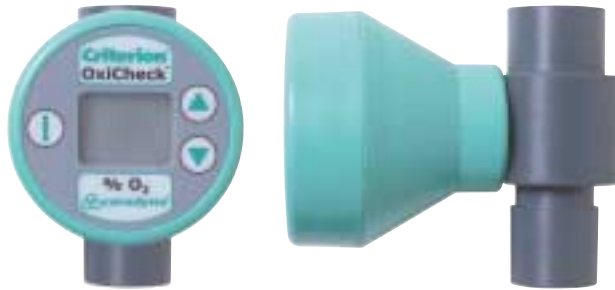

Criterion OxiCheck
Portable O₂ Analyzer

Product Specifications

Dimensions: 80x80x56mm

Weight: 90g

Performance Specifications

Measurement Range: 0-100%

Measurement Resolution: 1%

Analyzer Life: 500,000 oxygen hours
(15 months continuous use at 45% O₂, or 5000 hours at 100% O₂)

Operating Temperature: 0-50°C

Operating Humidity: 0-99% RH

Storage Temperature: -20° to 60°C

Accuracy: +/- 3%

Product Description

Part Number:	Type:	Qty/Box:
8-100661-00	OxiCheck Analyzer	1

The Criterion® OxiCheck™ is a compact analyzer designed to “spot check” oxygen concentration.

Designed to be quickly and easily connected and removed directly from a 22mm output connector of a Whisperflow® Generator or most ventilators, OxiCheck can be used to confirm the FiO₂ output being delivered to a patient very accurately.

It is small enough to be carried in a pocket and used with multiple patients.



Caution: Federal Law (USA) restricts this device to sale by or on the order of a licensed physician.



Caradyne is a registered trademark.
Criterion is a registered trademark and OxiCheck is a trademark of Caradyne, Ltd.



Manufactured by:
Caradyne
Parkmore Business Centre
Parkmore West, Galway, Ireland
T: + 353 (0) 91 709010
F: + 353 (0) 91 758929
Email: info@caradyne.com
Web: www.caradyne.com

Caradyne
2055 Executive Drive
Indianapolis, IN 46241 USA
T: 317.246.9010
F: 317.246.9011

Lit: xxxx Rev. A, 06/21/02